




Repairing nose pads on eyeglasses

Fix a loosen nose pad by re-screwing it back on the pad arm.

Written By: Vivi Chou

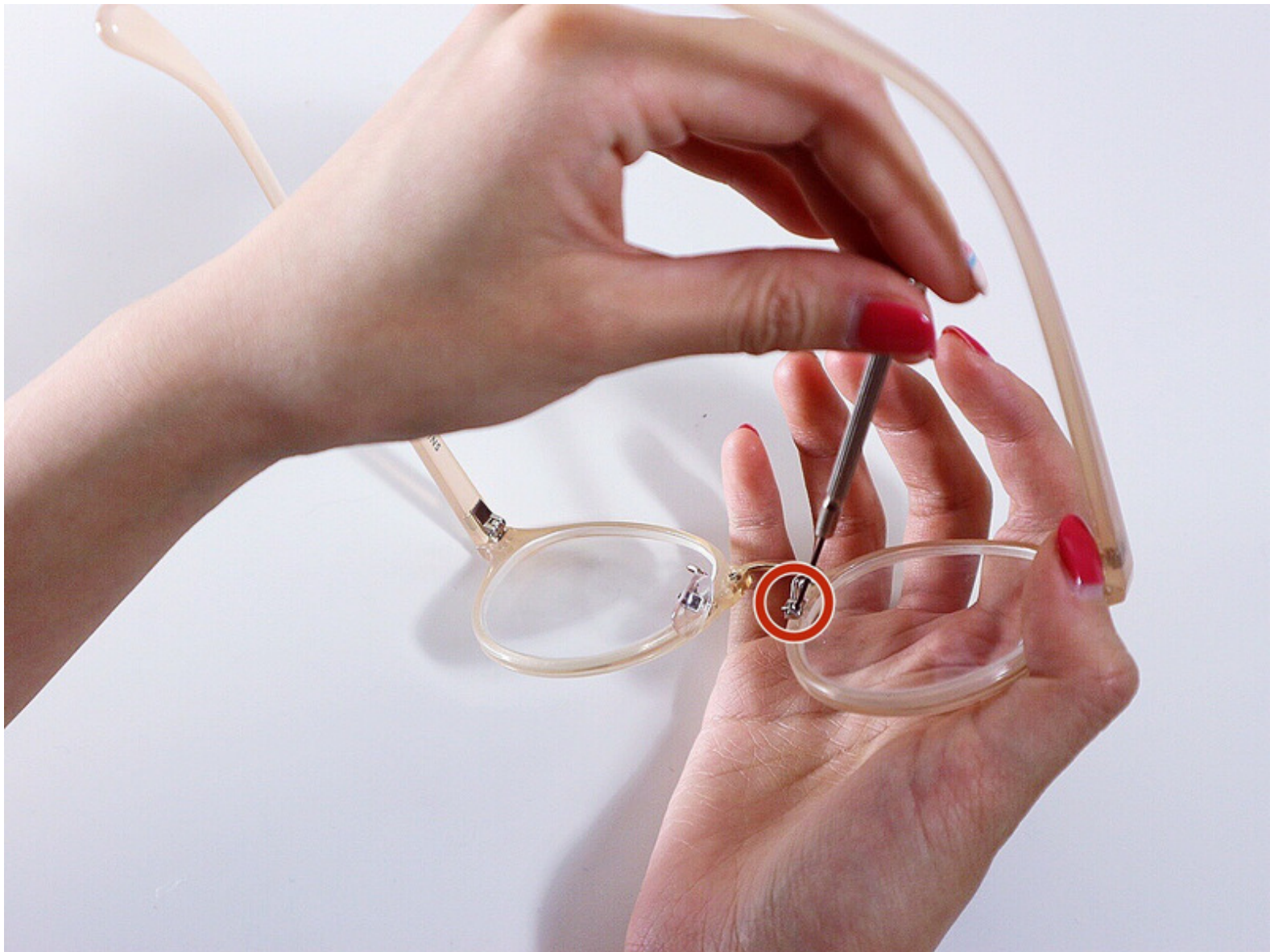

INTRODUCTION

This guide will walk you through our repair process step by step on how to reattach a loosen nose pad on a pair of eyeglasses. You will need a nose pad, a 1.0mm flathead screwdriver, a magnet, and an extra screw if you have lost the one on your eyeglasses.

We feel this guide is necessary because according to our research, myopia is affecting nearly 40 percent of the US population. In addition, based on our own experiences, nose pads are normally the weakest part of the glasses. For those who need their glasses everyday, knowing how to repair the nose pads by themselves is extremely useful.



TOOLS:



- [1.0 mm Flathead Screwdriver](#) (1)
- [Magnetic Telescoping Pick-Up Tool](#) (1)
- [Eyeglasses Repair Kit](#) (1)



PARTS:

- [Eyeglass Nose Pads](#) (1)
-


Step 1 — Repairing nose pads on eyeglasses

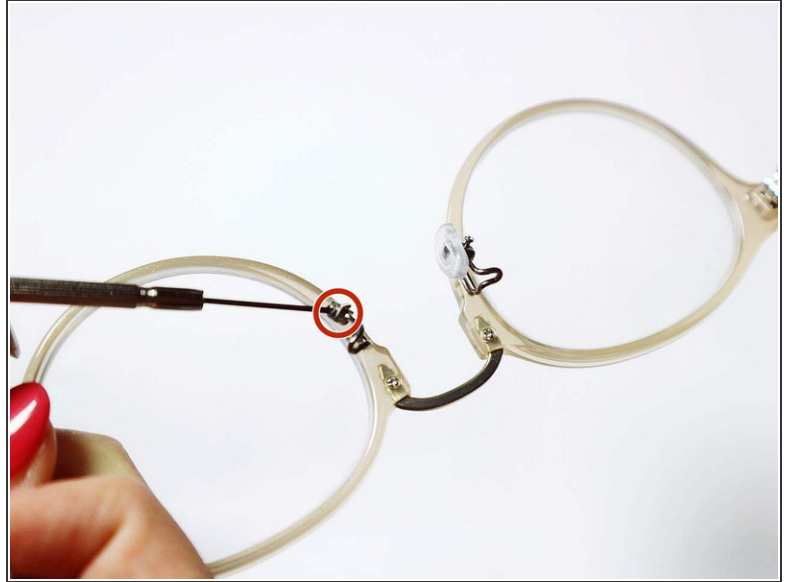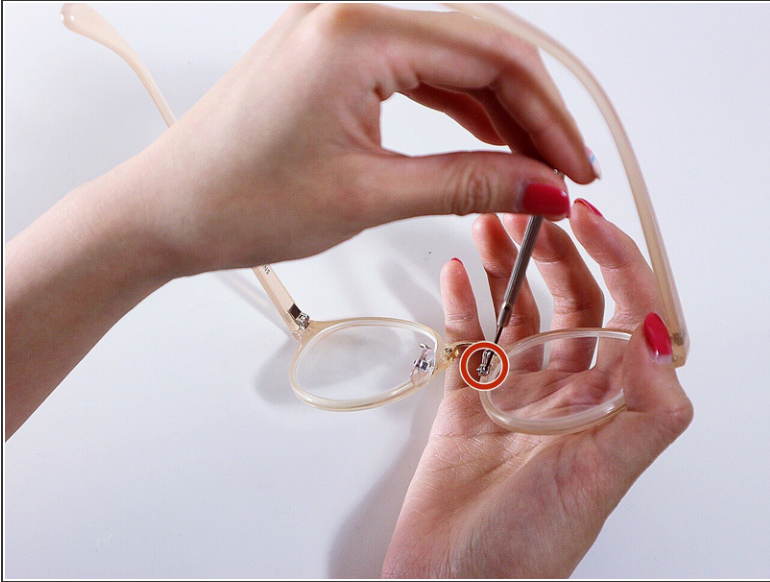

- Remove the screw from the pad arm.

ⓘ If the screw is missing, get a backup screw, and proceed to step 2.

Step 2

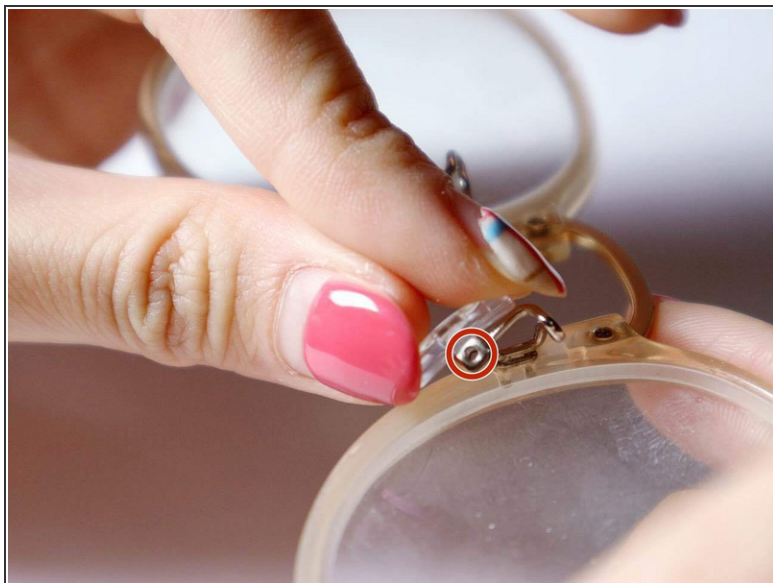



- Insert the nose pad into the mount on the frame.
- ⓘ Make sure the nose pad lines up with the screw hole.
- While holding the nose pad in place, use the magnet to attract the screw head.
- ⓘ Make sure the screw is vertical to the magnet and the screw head is above.

Step 3




- Place the screw into the screw hole.

Step 4

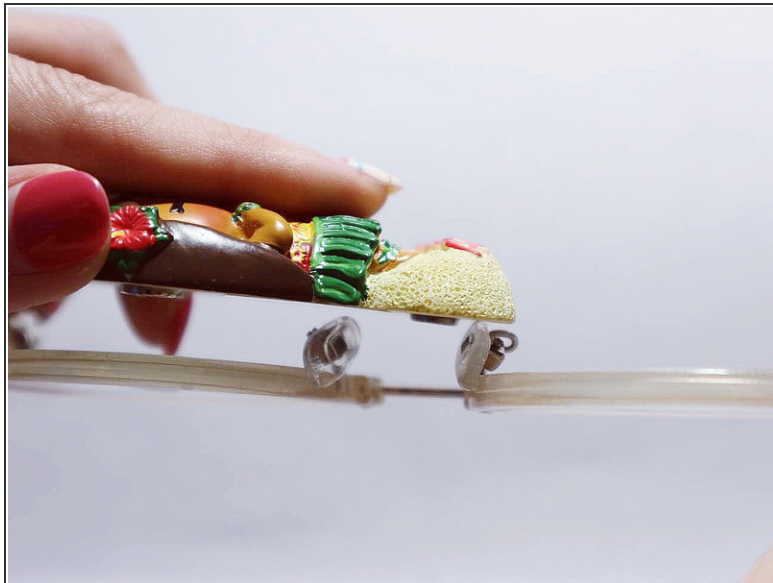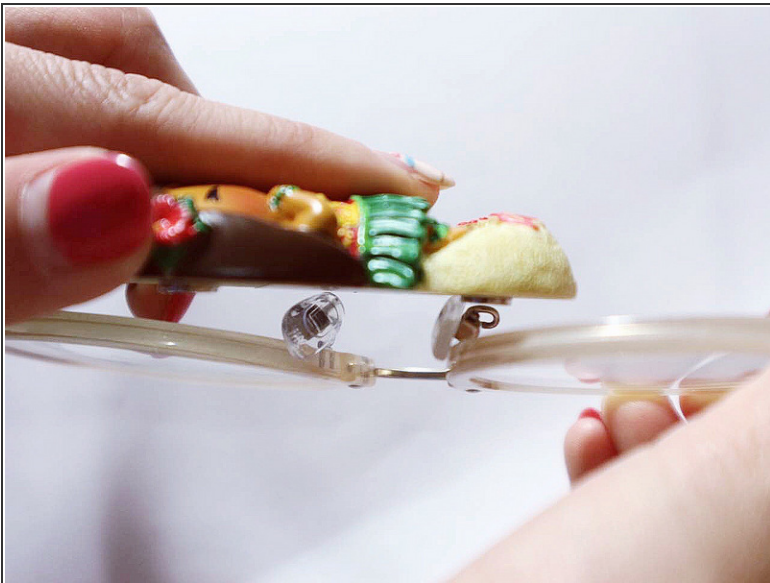

- Using a horizontal motion, remove the magnet.
- ⓘ Don't use a vertical motion to remove the magnet, otherwise the screw will come out.

Step 5

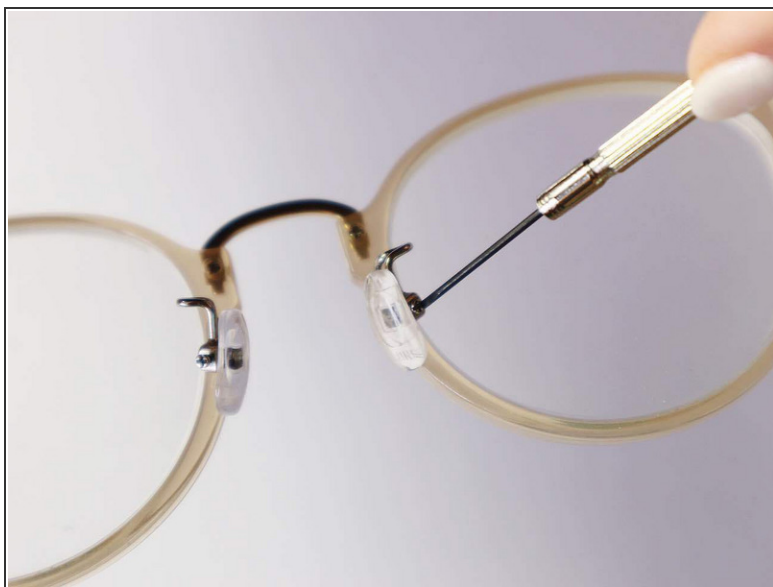

- Holding the glasses with the screw head pointing up, tight the screw with the screwdriver.
- ☑ Make sure the screw is tightened enough.

To reassemble your device, follow these instructions in reverse order.

This document was last generated on 2017-07-11 02:47:47 AM.